



# Diabetes Care




# Insulin Syringe

Syringes for subcutaneous insulin therapy injection.

## Feature:

- Premium needle quality and silicon-coated to ease skin penetration and minimize pain.
- High transparency barrel with clearly readable graduated scale for better dose precision.
- No-dead-space.

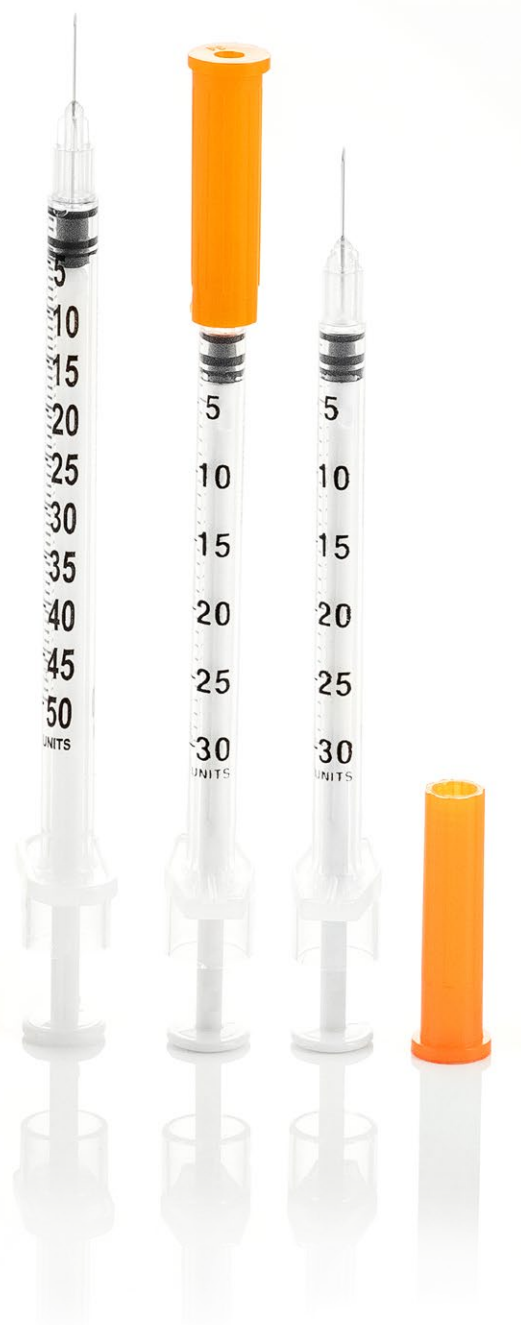

Code	Unit	Size	Needle size	Cap color
IS03	U-100	0.3ML	25-33GX5/16"-1/2"	orange
IS05	U-40/ U-100	0.5ML	25-33GX5/16"-1/2"	U-40 red/ U-100 orange
IS10	U-40/ U-100	1ML	25-33GX5/16"-1/2"	U-40 red/ U-100 orange
Needle sizes according to customer's requests				

# Safety Insulin Syringe

## Feature:


- Safety insulin syringes with integrated needle
- The innovative safety insulin syringe offers comfort and easy insulin therapy injection.
- Help eliminate accidental needle stick injury.



Diabetes Care

## Type 1

- Pull back the plunger after injection, the locking ring secures the needle in the barrel of the syringe protecting from needle stick injury.

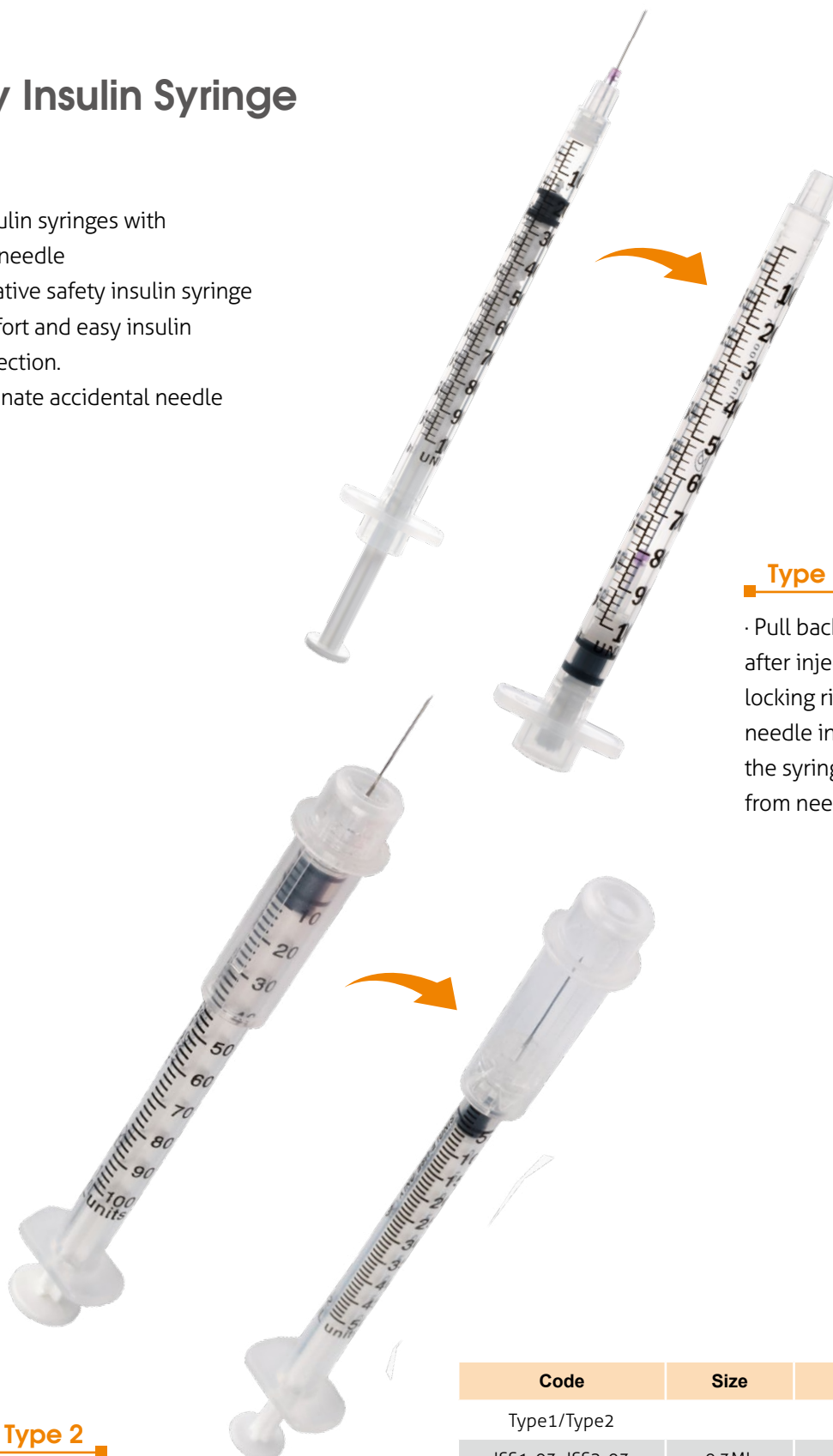

## Type 2

- Pull the safety cartridge to withdraw the needle into cartridge after injection.
- Easy to lock by rotating the safety cartridge.

Code	Size	Needle size
Type1/Type2		
ISS1-03 ISS2-03	0.3ML	29-33G*5/16"-1/2"
ISS1-05 ISS2-05	0.5ML	29-33G*5/16"-1/2"
ISS1-10 ISS2-10	1ML	28-33G*5/16"-1/2"



## Insulin Pen Needle

### Feature:

- Extremely fine three-facet grinding minimizes tissue trauma
- Special silicone finish for smooth gliding and virtually pain-free penetration of skin and tissue
- Extra-thin-wall technology allows for a larger inner diameter and an optimized Insulin flow; less force is required and injection time is reduced
- Compatible with all major pen devices available on the market



## Safety Insulin Pen Needle

### Feature:

- Automatically achieve the function of protecting needle tip after use.
- Easy self-disable set
- Easy to see the droplets
- Compatible with all pen injectors



Code	Gauze	size and length (mm)
Safety/Non-safety		
IPNS33 IPN33 <b>new</b>	33G	0.20X4-12mm
IPNS32 IPN32	32G	0.23X4-12mm
IPNS31 IPN31	31G	0.25X4-12mm
IPNS30 IPN30	30G	0.3X4-12mm
IPNS29 IPN29	29G	0.33x4-12mm

## Product Features



### Safety design and easy for using

The strength design in needle set make it more convenient and safe to use

### High compatibility and cover almost all branches

The product has high compatibility. The standard design of internal threaded ensure it can stably match all branches insulin pen.

### Best penetration make injection more comfortable


Insulin pen needle has a series of advanced and mature needle tip technique and silication technique. The sharp needle tip make puncture more easily and let patients feel more comfortable during injection.



Diabetes Care



## 33G, a shorter and finer insulin pen needle!!

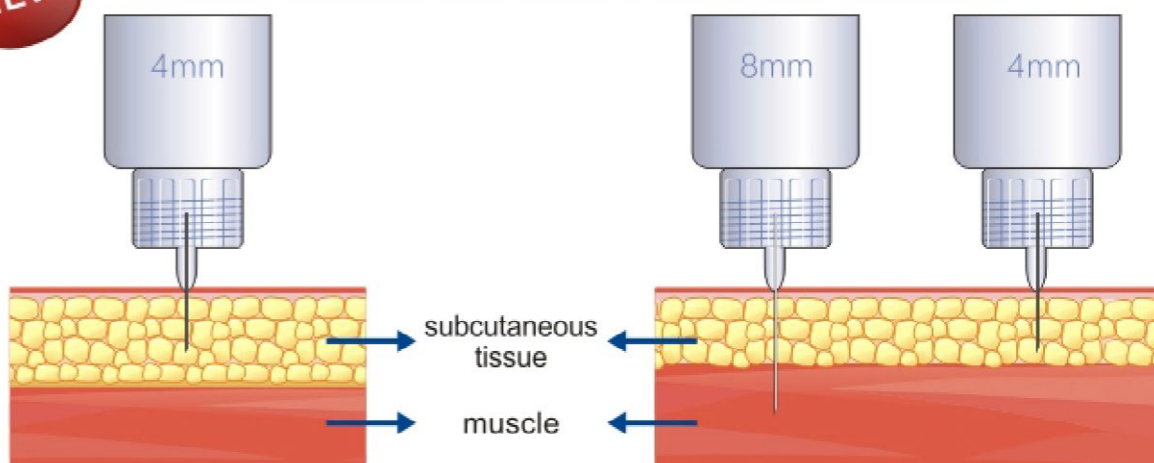

Don't need to pinch the skin when use 4mm needle

Contrast of 8mm, 4mm needle used in few subcutaneous tissue area

## High Matching Degree



NovoPen® 4



NovoPen® 3



NovoPen® Junior



NovoRapid® Flexpen®



BD: Puregon Pen®



Sanofi: Aventis Lantus



Sanofi: OptiPen®



Sanofi: TactiPen®



Owen Mumford: AutoPen®



Owen Mumford: AutoPen®2



Lily: HumaPen®



Lily: HumaPen®



Wanbang: Wangbang Pen 2



Tonghua dong bao: Doing Bao Pen®



Tonghua dong bao: Gansulin Pen®

- The screw thread design of insulin pen needle meets the standard regulated by ISO11608-2 which is compatible with other standard ones domestically produced or imported.



## Safety Lancet

It is used in blood collection sampling and help eliminate accidental needlestick injuries.

### Feature:

- Needle is fully shielded before and after use
- Cannot be reused
- Ultra sharp tri-bevel needle, high speed penetration
- Multiple depth selection with color coding



Type3



Type2

Type1



Type5



Type4

Code					Size	Length
TYPE 1	TYPE 2	TYPE 3	TYPE 4	TYPE 5		
SBL01-30	SBL02-30	SBL03-30	SBL04-30	SBL05-30	30G	1.8mm
SBL01-28	SBL02-28	SBL03-28	SBL04-28	SBL05-28	28G	1.8mm
SBL01-26	SBL02-26	SBL03-26	SBL04-26	SBL05-26	26G	1.8mm
SBL01-23	SBL02-23	SBL03-23	SBL04-23	SBL05-23	23G	1.8mm
SBL01-21	SBL02-21	SBL03-21	SBL04-21	SBL05-21	21G	1.8/2.2/2.8mm
SBL01-18	SBL02-18	SBL03-18	SBL04-18	SBL05-18	18G	1.8mm



# Blood Glucose Monitoring Systems

## Feature :

- Fast, Safe and Convenient.
- Small blood sample requirement reduces pain & sensitivity for easy diabetes monitoring.
- Self-testing on the volume of the blood sample
- Coding/no-coding optional



Small blood sample  
Less pain



Alarm



5 seconds  
Fast testing



Auto-off



Strip ejector



High accuracy  
ISO 15197:2013  
certified



No coding



Test strip shelf  
life 2 yearsTest



360 test memories



## BG205

### Feature :

- 1ul small blood sample
- 5seconds measuring time
- 360 reading memories
- Extra large LCD
- Backlight optional
- Strip ejector
- 7,14,28 day average
- Alarm clock reminder
- Switchable dual scale
- Auto Power off
- Deluxe carry case
- Coding/no-coding optional



Diabetes Care



## BG202

### Feature :

- 1ul small blood sample
- 5 seconds measuring time
- 360 reading memories
- 7,14,28 day average
- Alarm clock reminder
- Switchable dual scale
- Auto Power off
- Deluxe carry case
- Coding/ no-coding optional



## BG203

### Feature :

- 1ul small blood sample
- 5 seconds measuring time
- 360 reading memories
- Extra large LCD
- Backlight optional
- Strip ejector
- 7,14,28 day average
- Alarm clock reminder
- Switchable dual scale
- Auto Power off
- Deluxe carry case
- Coding/ no-coding optional

## BG208

### Feature :

- 1ul small blood sample
- 5 seconds measuring time
- Strip ejector
- 360 reading memories
- 7,14,28 day average
- Alarm clock reminder
- Switchable dual scale
- Auto Power off
- Deluxe carry case
- Coding/nocoding optional





## Accessories



**Sharp disposal container**

0.2L/1Quart/1L/2Quart/ 2L



**Alcohol swab**

3x6cm/ 4x8cm



**Lancet**

23g /26g /28g /30g /32g /33g



**Lancing devices**

The adjustable tip offer 5 depth levels



**Safety Blood Lancet**

18g/21g/23g /26g /28g /30g



**Test strips**

Width 5.5mm/6.5mm

